
Oxygenator with integrated arterial filter / Oxygénateur avec filtre  
artériel intégré / Oxygenator mit integriertem arteriellem Filter  
Oxigenador con filtro arterial integrado / Ossigenatore con filtro  
arterioso integrato / Oxygenator met geïntegreerde arteriële filter

**CAPIOX<sup>®</sup> FX 05**  **BABY-FX<sup>™</sup>**

*Applicable for Vacuum Assisted Venous Drainage  
Utilisable pour Drainage Veineux Assisté sous Vide  
Einsetzbar für Vakuumunterstützte Venen-Drainage  
Aplicable al Drenaje Venoso Asistido por Vacío  
Utilizzabile per il Drenaggio Venoso Assistito con il Vuoto  
Toepasbaar voor Veneuze Drainage Onder Vacuüm*

**Instructions for Use / Mode d'emploi  
Bedienungsanleitung / Instrucciones de Uso  
Istruzioni per l'Uso / Gebruiksaanwijzing**

**Xcoating<sup>™</sup>** **CE** 0197

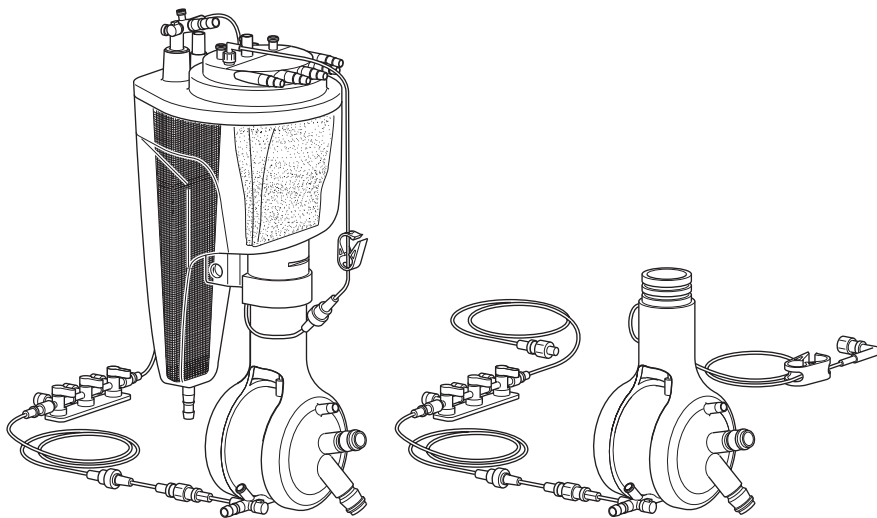**REF : CX\*FX05RW  
CX\*FX05RE**

**REF : CX\*FX05W  
CX\*FX05E**

**TERUMO<sup>®</sup>**

ENGLISH

FRANÇAIS

DEUTSCH

ESPAÑOL

ITALIANO

NEDERLANDS

FX35M051


**SYMBOLS ON THE PACKAGES OF THIS PRODUCT / SYMBOLES SUR LE  
CONDITIONNEMENT DU PRODUIT / AUF DER VERPACKUNG DES PRODUKTES  
ABGEBILDETE SYMBOLE / SIMBOLOS UTILIZADOS EN EL ENVASE DE ESTE PRODUCTO /  
SIMBOLI PRESENTI SUL CONFEZIONAMENTO DI QUESTO PRODOTTO / SYMBOLEN  
AANWEZIG OP DE VERPAKKING VAN DIT PRODUCT**


Code No. Référence Artikelnummer Código Codice Artikelnummer	Lot No. Numéro de lot Chargen-Bezeichnung Lote Lotto Lot	Expiry date Utiliser avant Verwendbar bis Caducidad Scadenza Vervaldatum	Sterilized by ethylene oxide. Stérilisé à l'oxyde d'éthylène. Sterilisiert mit Ethylenoxid. Esterilizado con óxido de etileno. Sterilizzato con ossido di etilene. Gesteriliseerd met ethyleenoxide.	Do not reuse. Ne pas réutiliser. Nicht zur Wiederverwendung. No reutilizare. Non riutilizzare. Niet opnieuw gebruiken.	Read the instructions for use. Lire le mode d'emploi. Lesen Sie die Gebrauchsanweisungen. Leer las instrucciones de uso. Leggere le istruzioni per l'uso. Lees de gebruiksaanwijzing.
---	---	---	---	---	--

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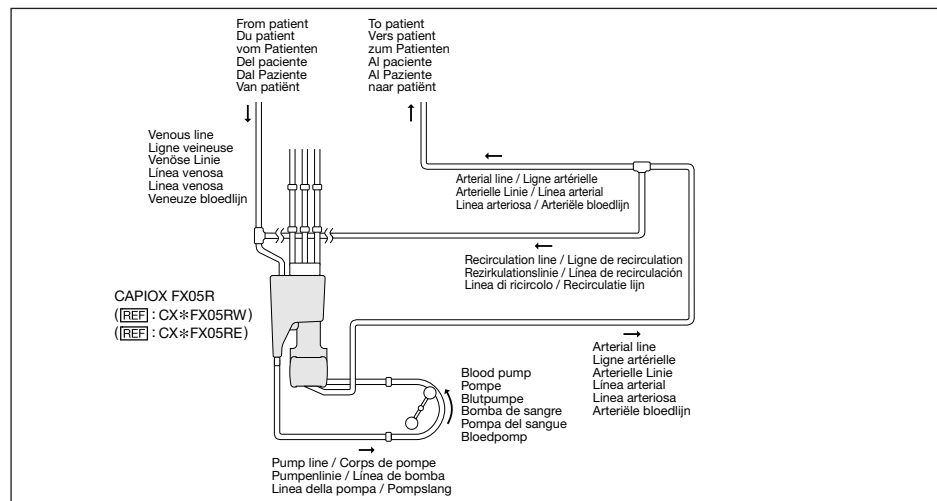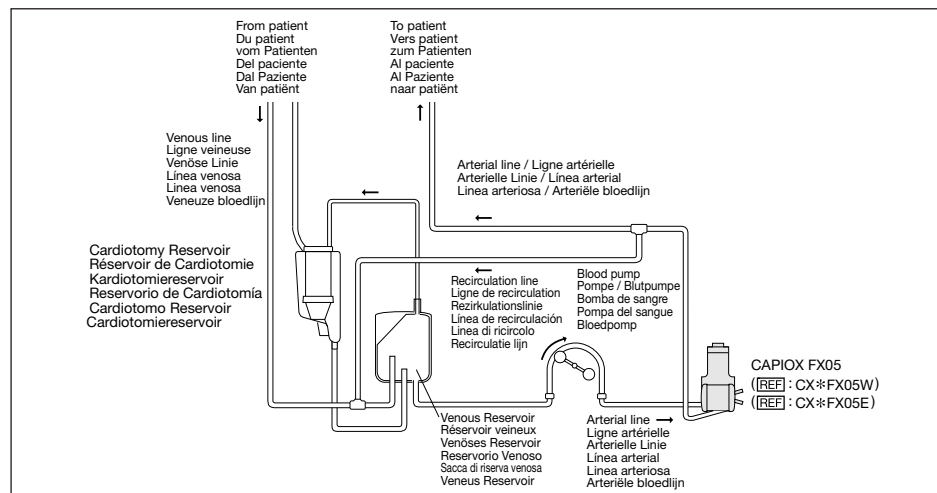**CAPIOX® FX 05**  
**BABY-FX™**

**Applicable for Vacuum Assisted Venous Drainage**

**Instructions for Use**  
**P.3-P.16**

Read all warnings, cautions and instructions carefully prior to use.

**DESCRIPTION OF PRODUCTS**

The CAPIOX FX05 is a membrane oxygenator with microporous polypropylene hollow fibers, to be used as an extracorporeal gas exchange device in which the blood flows on the outside of the fibers and with the ventilating gas flowing through the fibers.

The CAPIOX FX05 consists of a module for gas exchange with an integrated heat exchanger and Arterial Filter. The CAPIOX FX05 is available with a Hardshell Reservoir with an integrated cardiomy filter to provide a single-unit construction for ease of use (Code No.: CX \* FX05RW, CX \* FX05RE) or without a Hardshell Reservoir (Code No.: CX \* FX05W, CX \* FX05E).

(CX \* FX05RW, CX \* FX05W and CX \* FX05RE, CX \* FX05E differ in the positioning of the arterial blood port of the oxygenator.)

The blood contacting surfaces are coated with **Xcoating**.

**Xcoating** is a bio-compatible material that is applied on the blood-contacting surfaces of the device to reduce the adhesion of platelets to the device.

**INDICATION FOR USE**

The CAPIOX FX05 is intended to be used during open heart surgical procedures to transfer oxygen and remove carbon dioxide from blood and to control the blood temperature during cardiopulmonary bypass for periods up to 6 hours. The CAPIOX FX05 is a NEONATE, INFANT oxygenator intended for use in procedures up to maximum flow of 1.5 L/min. The patient weight and BSA should be considered upon use.

The FX05 Hardshell Reservoir is also intended for use in vacuum assisted venous drainage procedures. (Refer to "Vacuum Assisted Venous Drainage", page 12)


The Integral Arterial Filter is intended to filtrate non-biologic particles and emboli and to facilitate gaseous emboli removal from the blood flowing through a cardiopulmonary bypass circuit.

**DIRECTION OF BLOOD OUTLET PORT**


The oxygenator is available with TWO OUTLET PORT CONFIGURATIONS - LEFT side (W) or RIGHT side (E). Refer to the diagram below to determine which configuration is needed for your circuit.

TYPE	LEFT	RIGHT
Oxygenator with Hardshell Reservoir	code : CX*FX05R <b>W</b>	code : CX*FX05R <b>E</b>
Oxygenator without Hardshell Reservoir	code : CX*FX05 <b>W</b>	code : CX*FX05 <b>E</b>
Place waterports facing South as indicated in the illustration. The blood outlet port orientation is (W)est - LEFT side or (E)ast - RIGHT side.		


## STRUCTURE

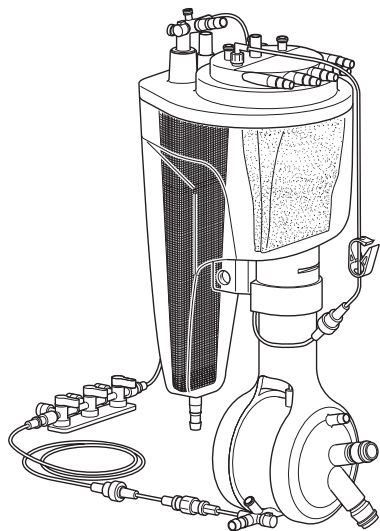
REF : CX\*FX05RW

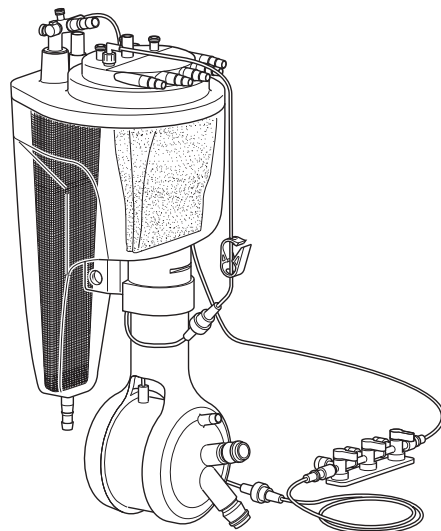REF : CX\*FX05RE

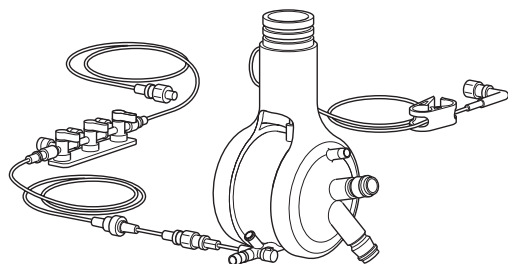
REF : CX\*FX05W

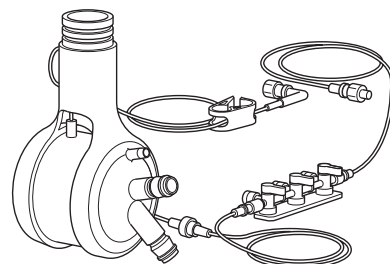
REF : CX\*FX05E

CX \* FX05RW, CX \* FX05W and CX \* FX05RE, CX \* FX05E differ in the positioning of the arterial blood outlet port of the oxygenator.  
All figures hereafter only refer to FX05RW.

## STRUCTURE

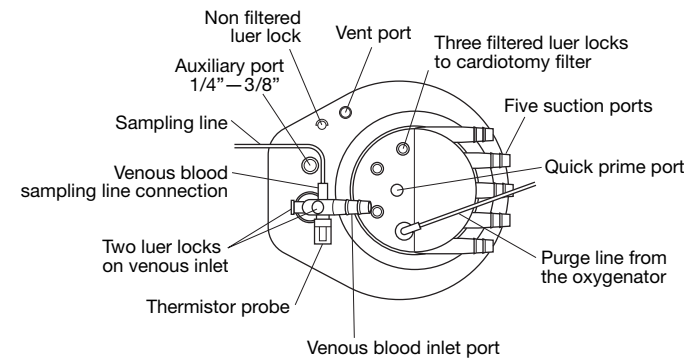Hardshell Reservoir
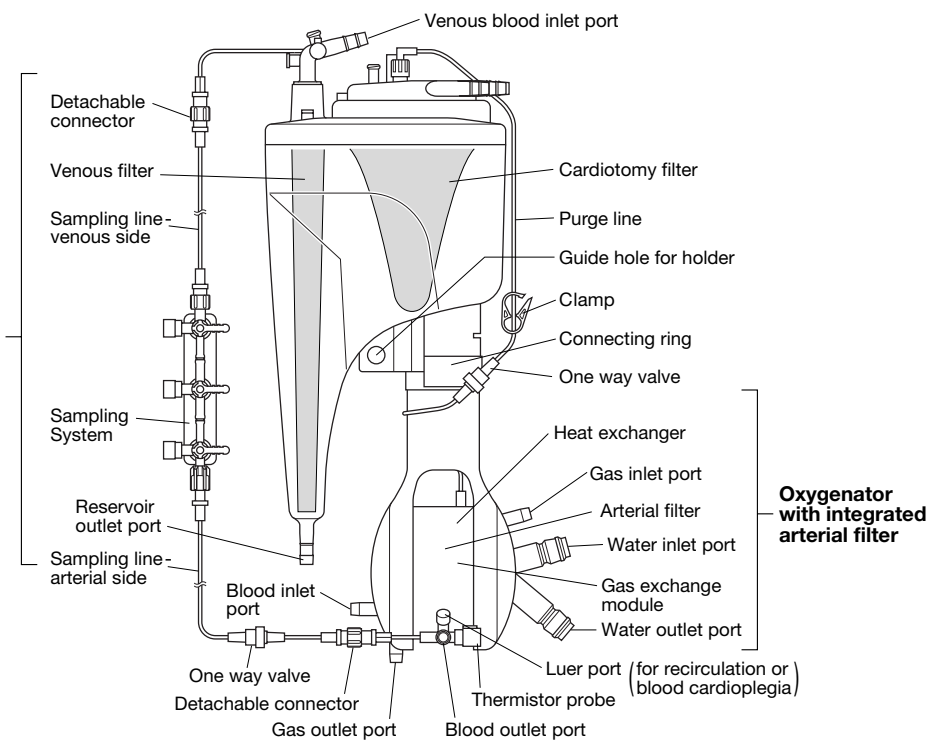
## WARNINGS

Describe serious adverse reactions and potential safety hazards, limitations in use imposed by them, and steps that should be taken if they occur.

### WARNINGS

- CAPIOX FX05 is designed to operate at blood flow rates within the range of 0.1 to 1.5 L/min. Do not use any blood flow rate outside this range.
- Do not use solvents such as alcohol, ether, acetone, etc. These solvents may cause damage if used in or on the device.
- In order to prevent gaseous emboli from entering the blood phase, follow the instructions below.
  - DO NOT OBSTRUCT GAS OUTLET PORT.
  - Pressure in the blood phase should always be higher than that in the gas phase.
  - The gas flow rate should not exceed 5 L/min. Excessive gas flow rate will bring about pressure increase in the gas phase.
  - During recirculation, do not use pulsatile flow and do not stop the blood pump suddenly as these actions may cause gaseous emboli to enter the blood phase from the gas phase due to inertia force.
  - When CAPIOX FX05 Oxygenator module is used separately from the Hardshell Reservoir, set the module so that the upper end of the fibers is lower than the blood level in the venous reservoir. This prevents gaseous emboli from entering the blood phase from the gas phase.
  - To prevent gaseous emboli from entering the blood phase, make sure that the arterial pump flow rate always exceeds the flow rate of the cardioplegia line. The blood flow rate of the cardioplegia line should not exceed 0.5 L/min.
  - Total flow rate of the arterial line and any separate arterial lines must not exceed the flow rate at the oxygenator inlet port.
  - Do not exceed a 15°C (27°F) temperature difference between blood and water in heat exchanger to prevent gas dissolved in blood from forming bubbles.
- Pressure at the blood inlet of the oxygenating module should not exceed 133 kPa (1,000 mmHg). Pressures greater than 133 kPa (1,000 mmHg) may cause leaks or damage to the device.
- Water pressure at the heat exchanger inlet should not exceed 196 kPa (2 kgf/cm<sup>2</sup>). Pressures greater than 196 kPa (2 kgf/cm<sup>2</sup>) may cause leaks or damage to the device.
- Adequate heparinization of the blood is required considering patient condition and perfusion technique to prevent it from clotting in the system.
- Stop the gas flow when circulation is suspended. During recirculation, check blood gas pressure. Excessive gas flow may cause low PaCO<sub>2</sub>, alkalosis, or blood damage.
- Use caution when removing air during priming and perfusion. Excessive shock to the device - especially with hard objects - can cause damage to the device.

Use CAPIOX FX05 Oxygenator with Hardshell Reservoir taking note of the additional warnings below.

- The blood flow into the cardiotomy filter should not exceed the rate of 1.5 L/min. Excessive blood flow rate may increase the pressure in the cardiotomy filter, resulting in back-flow into any solution or blood administration line connected to the Hardshell Reservoir.
- Yellow vent port cap need not be removed, since this port assures adequate venting with the cap attached. Do not close the vent port, as this may cause positive pressure in the Hardshell Reservoir, resulting in back-flow into solution or blood administration lines connected to the Hardshell Reservoir.
- Negative pressure below -20 kPa (-150 mmHg) should not be applied to this reservoir as this may damage it.
- Minimum operating volume in the reservoir is 15 mL. Volume below 15 mL will pull air from the reservoir. Set appropriate blood storage level, relative to venous flow rate, to prevent gaseous emboli passing to patient.  
(Refer to Fig.10, page 14 "Safety Margin")

## PRECAUTIONS

Include information regarding any special care to be exercised by the practitioner for the safe and effective use of the device.

### PRECAUTIONS

- This device should only be used by properly trained and qualified personnel.
- This device is sterile and non-pyrogenic in the unopened, undamaged package. Inspect the device and package carefully. Do not use if the package and / or device is damaged, or if caps are not in place.
- Dispose of safely after single use to avoid risk of infection.
- DO NOT RESTERILIZE.
- This device is intended for single use only. DO NOT REUSE as reprocessing may compromise the biocompatibility and functional integrity of this product.
- In case of the need to replace the oxygenator, make sure that a spare one is available.
- Use aseptic technique in all procedures.
- Perform priming using crystalloid solution which contains no blood, plasma and/or blood derived products. If blood derived products are used during priming, start bypass after performing air purge thoroughly.
- This device should be used with an appropriate holder supplied by TERUMO.
- Place CAPIOX FX05 Oxygenator module downstream from the blood pump.
- Thermistor probes can be used with a temperature monitor from Measurement Specialties Inc. by use of the cables supplied.
- Temperature of water delivered to heat exchanger should not exceed 42°C (108°F), otherwise blood damage may result.
- A minimum of 0.05 L/min oxygen gas flow or minimum V/Q 0.2 is needed when blood is circulated. Less than 0.05 L/min oxygen gas flow or V/Q less than 0.2 may result in inadequate gas exchange.
- Use of a pre-bypass filter is recommended to trap any particulate matter in the circuitry and priming solution.
- The Sampling Line is detachable. If the line is detached and not in use, make sure to cap the end.
- Use of a recirculation line is recommended in the arterial line after the FX oxygenator to facilitate gaseous emboli removal during priming and emergency situations.
- It is recommended to use an air bubble detector between the reservoir outlet and FX oxygenator inlet. In the event gaseous emboli are detected, determine the source and take corrective action to eliminate further introduction of gaseous emboli into the perfusion circuit.
- Ensure the oxygenator is never subjected to a negative pressure. Negative pressure in the oxygenator can cause gaseous emboli to enter the blood phase from the gas phase.

Use CAPIOX FX05 Oxygenator with Hardshell Reservoir taking note of the additional precautions below.

- When rotating the oxygenator module, be sure to hold the lower part of the connecting ring. If the oxygenator module is rotated while holding the oxygenator body, it could be damaged.
- When rotating the venous blood inlet port on the header of Hardshell Reservoir, make sure that venous line and thermistor probe cable do not come into contact with any solution or administration line connected to the Hardshell Reservoir to prevent it from kinking.
- When rotating the suction port on the header of Hardshell Reservoir, make sure that the purge line does not come into contact with any solution or administration line connected to the Hardshell Reservoir to prevent it from kinking.
- When separating Hardshell Reservoir from oxygenator module, hold the oxygenator firmly and remove connecting ring.
- The caps of unused ports should be left in place. This avoids contamination and prevents leakage of blood.
- Make sure the caps of unused luer ports are tightened firmly to prevent the leakage.
- Cardiotomy filter must be wetted prior to use to achieve maximum flow capacity.
- Where appropriate, introduce blood, plasma and/or blood derived products through quick prime port or any of the luer ports leading to cardiotomy filter after debubbling of oxygenator.
- Scientific literature contains reports of altered dose response to some medications, such as nitroglycerin, fentanyl, etc., possibly due to varying degrees of absorption by synthetic materials such as those found in extracorporeal circuits.
- For a low dose of drug administration, do not inject from luer ports leading to the cardiotomy filter. If drug is injected inside the cardiotomy filter, it could stagnate in the filter.
- When storing the product, use caution not to get the product wet, and avoid any exposure to direct sunlight or high humidity.

CAUTION : Federal law (U.S.A.) restricts this device to sale by or on the order of a physician.

## METHOD OF OPERATION

Read WARNINGS (page 6) and PRECAUTIONS (page 7) before use.  
The following is a basic description of the use of the CAPIOX FX05 Oxygenator and Hardshell Reservoir.

### a. SET-UP

1. Remove CAPIOX FX05 from package and examine it for defects.

**CAUTION** Do not use if the package or device is damaged (e.g. cracked) or any of the port caps are off.

2. Place CAPIOX FX05 in its holder.

- a. CAPIOX FX05 Oxygenator with Hardshell Reservoir:

Set the holder at the appropriate height to prevent the arterial line and other lines from kinking.

Making sure that the connecting ring is in place, set oxygenator in its holder (Code No.:XX \*CXH05R) with lever out, then push the lever in. (See Fig.1) If the oxygenator is not set into the holder correctly, the amount of fluid in the reservoir cannot be correctly assessed.



- b. CAPIOX FX05 Oxygenator module without Hardshell Reservoir:

Set the holder at the appropriate height to prevent the arterial line and other line from kinking.

Set oxygenator module in its holder (Code No.:XX \*CXH05) with lock cover raised, then close the cover. (See Fig.2)

NOTE : • When using CAPIOX FX05 Oxygenator module and Hardshell Reservoir separately,

- a) Hold the module firmly and separate it from Hardshell Reservoir by removing connecting ring.
- b) Set reservoir in the holder (Code No.:XX \*CXH05R) for CAPIOX FX05 Oxygenator with Hardshell Reservoir and push the lever in.
- c) Set the oxygenator module in its holder (Code No.:XX \*CXH05) and close the "lock" cover.

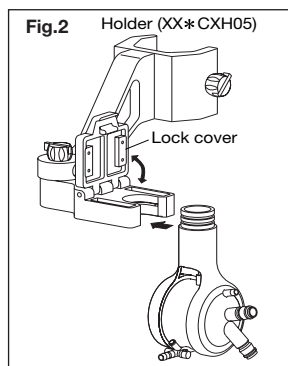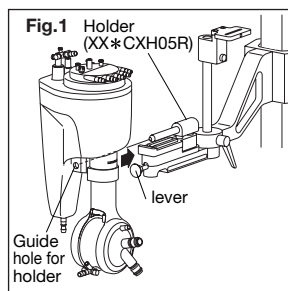**CAUTION** If the product is dropped during set-up, do not use it. Replace with another device.

3. Connect water lines (1/2"(12.7mm) tubing or couplers) to the water ports of CAPIOX FX05. (See Fig.3 - (1) )

**CAUTION** Use the upper port for water supply and the lower port for water drainage. Otherwise the heat exchange will not be performed sufficiently.

4. Start water circulation through heat exchanger and circulate for at least 5 minutes. Check for leaks.


**WARNING** Do not use an oxygenator that leaks.

5. Attach 1/4"(6.4mm) venous line to the venous blood inlet port of Hardshell Reservoir. (See Fig.3 - (2) ) If 3/8"(9.5mm) venous line is used, connect the included 3/8" adapter to the venous blood inlet port and connect the 3/8" adapter to the venous line.

6. Attach 1/4"(6.4mm) arterial line to the blood outlet port of CAPIOX FX05 Oxygenator. (See Fig.3 - (3) ) If 3/16"(4.8mm) arterial line is used, connect the included 3/16" adapter to the arterial blood outlet port and connect the 3/16" adapter to the arterial line.

7. Remove the cap from the luer port on the oxygenator outlet port (See Fig. 3 - (4) ) and attach the recirculation line. Connect the other end to luer port on reservoir.

Fig.3

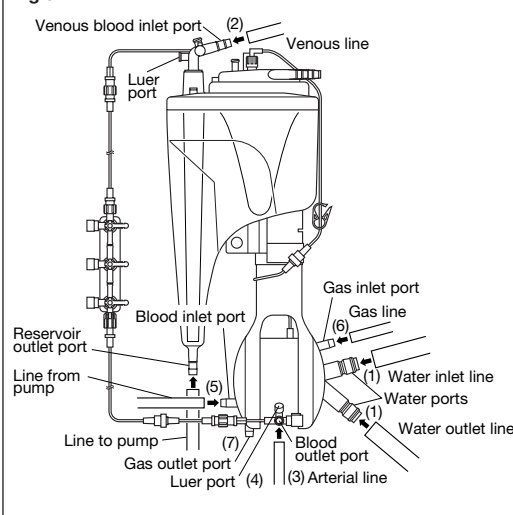**WARNING** If the port is not to be used, push cap in and twist it to seal tightly.

8. Connect 1/4"(6.4mm) pump line to the reservoir outlet port. Connect the other end to 1/4"(6.4mm) blood inlet port of oxygenator. (See Fig.3 - (5) ) If 3/16"(4.8mm) pump line is used, attach the included 3/16" adapters to the reservoir outlet port and to the blood inlet port of oxygenator and connect the pump line to the adapters.

9. Band all connections in the circuit.

10. Connect 1/4"(6.4mm) gas line to gas inlet port. (See Fig.3 - (6) )

**WARNING** DO NOT OBSTRUCT GAS OUTLET PORT. (Fig.3 - (7))

11. Thermistor probes can be used with Y.S.I. 400 series\* temperature monitor by use of the following parts.

Blue Cable (venous line) : Code No.: CX \*BP021

Red Cable (arterial line) : Code No.: CX \*BP022

\* Y.S.I. 400 (trade name : Measurement Specialties Inc.)

12. Connect suction lines and vent line to suction ports on Hardshell Reservoir. Suction ports of CAPIOX FX05 Hardshell Reservoir have blue caps. Remove the blue caps before connection. (See Fig.4)

**CAUTION** The suction ports on the Cardiotomy Reservoir section step-down from 1/4"(6.4mm) to 3/16"(4.8mm). Insert 3/16"(4.8mm) tubing up to the second barb and 1/4"(6.4mm) tubing well-beyond the second barb to insure security of connections to these ports.

13. For use of CAPIOX FX05 Oxygenator module without Hardshell Reservoir, connect venous side male luer lock of sampling system to venous line. (See Fig.5)

14. When the luer port located on venous blood inlet port of CAPIOX FX05 Hardshell Reservoir is to be used, attach a three-way stopcock before use.

15. Check figures and labels to make sure tubing is connected to the correct port.

Fig.4


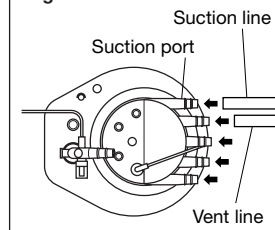
Fig.5
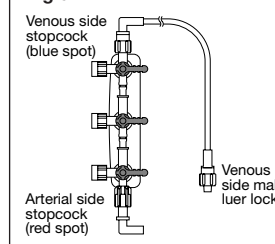
### b. PRIMING PROCEDURE

**CAUTION** Perform priming using crystalloid solution which contains no blood, plasma and/or blood derived products. Using blood derived products will increase the amount of time necessary to prime. Therefore, if blood derived products are used during priming, start bypass after performing air purge thoroughly.

NOTE : • CAPIOX FX05 can be primed without CO<sub>2</sub> flush when crystalloid priming solution is used. Use of CO<sub>2</sub> flush permits debubbling in a shorter time.

1. If recirculation line is used to prime, clamp arterial and venous lines proximal to recirculation line, and make sure that recirculation line is not clamped.
2. Introduce crystalloid priming solution through quick prime port or any of the luer ports leading to cardiotomy filter.
3. Make sure the circuit and the purge line are not clamped, then start pump at a low speed. After checking for leakage or any other problem, gradually increase flow above 0.5 L/min, but do not exceed 1.5 L/min. Vigorously recirculate the priming fluid through the entire circuit until all air bubbles are eliminated. After all air bubbles are eliminated, circulate at full flow for 10 min to check oxygenator and tubing for leakage or any other problem.

**WARNING** • Do not use an oxygenator and reservoir that leaks. Replace it with another CAPIOX FX05 oxygenator and reservoir.

• Do not use a tubing with internal diameter less than 3/16"(4.8mm) as a recirculation line. Also do not use a sampling line nor a purge line for recirculation. If used, the oxygenator module can be damaged as a result of excessive positive pressure being generated inside itself.

4. When appropriate after the debubbling, introduce blood or blood-derived products through quick prime port or any of the luer ports leading to cardiotomy filter.


- Set stopcocks as shown in Fig. 6, and close sampling line with arterial side stopcock to prevent arterial to venous shunting during extracorporeal circulation. After closing the purge line, reduce blood flow rate gradually to zero; then close recirculation line.

- After priming, if air bubbles continue to appear, identify the cause and make necessary corrections. Remove the air while opening the purge line.

**WARNING** During recirculation, do not use pulsatile flow or stop blood pump suddenly. Otherwise gaseous emboli may enter the blood phase from the gas phase due to inertia force.

**CAUTION** Close purge line before stopping bypass.
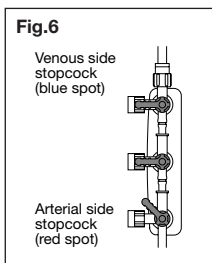
### c. INITIATION OF BYPASS

Initiate extracorporeal circulation using normal procedure, taking note of the following warnings.

- WARNINGS**
- Start gas supply only after blood circulation is initiated.
  - Before starting gas supply, confirm again that gas outlet port is not obstructed. Such obstruction may lead to pressure build-up in the gas phase, allowing gaseous emboli to enter into the blood phase.
  - Before initiating extracorporeal circulation, confirm that the recirculation line, purge line and sampling line are closed with arterial side stopcock. Otherwise, opening arterial line will cause the blood to flow back into reservoir through sampling line because of the patient's blood pressure and the head height.
  - Start gas supply with  $V/Q = 1$  and  $FiO_2 = 100\%$ , then make adjustments based on blood gas measurements.

### d. DURING PERFUSION

- In order to collect proper blood samples, withdraw at least 6 mL of blood, then collect blood through sampling line. In the case of arterial blood sampling, blood can be collected after opening the stopcock for arterial-venous shunting through sampling line.

**WARNING** Collect blood only while pump is running, or blood-side pressure will decrease and air bubbles may result.

NOTE : • For the use of the sampling system separated from Hardshell Reservoir, sampling manifold holder (Code No.: XX \*XH051) is available.

- Measure blood gases and make necessary adjustments as follows.
  - Control  $PaO_2$  by changing concentration of oxygen in ventilating gas using gas blender.
    - To decrease  $PaO_2$ , decrease  $FiO_2$ .
    - To increase  $PaO_2$ , increase  $FiO_2$ .
  - Control  $Paco_2$  by changing the total gas flow.
    - To decrease  $Paco_2$ , increase total gas flow.
    - To increase  $Paco_2$ , decrease total gas flow.

**WARNING** A phenomenon called wet lung may occur when water condensation occurs inside fibers of microporous membrane oxygenators with blood flowing exterior to the fibers. This may occur when oxygenators are used for a longer period of time. If water condensation and/or a decrease in  $PaO_2$  and/or an increase in  $Paco_2$  is noted during extended oxygenator use, briefly increasing the gas flow rate may improve the performance. Increase gas flow rate, to 5 L/min for 10 seconds. DO NOT repeat this flushing technique, even if oxygenator performance is not improved.

**CAUTION** A minimum of 0.05 L/min oxygen gas flow or minimum  $V/Q$  0.2 is needed when blood is circulated. Less than 0.05 L/min oxygen gas flow or  $V/Q$  less than 0.2 may result in inadequate gas exchange.

- The minimum reservoir level is 15 mL. However, adjust blood flow rate, making sure that the blood level in the reservoir is adequate, to prevent gaseous emboli passing to patient. (When using CAPIOX FX05 with Hardshell Reservoir, refer to Fig.10, page 14 "Safety Margin".)

- For purging, open the purge line.

**CAUTION** Close purge line at termination of bypass.

- Do not give any shocks to the product during perfusion.

### e. TERMINATION OF BYPASS

Terminate extracorporeal circulation using normal procedure, taking note of the following warnings :

- WARNINGS**
- Before reducing blood flow rate, confirm that purge line is closed and that sampling line is also closed with arterial side stopcock.
  - When stopping the pump, stop gas flow at once.
  - If circulation is to be resumed, recirculation at a low flow rate is recommended. During recirculation, check blood gas pressure. An excessive gas flow may cause low  $PaCO_2$ , alkalosis, or blood damage.

### f. OXYGENATOR REPLACEMENT




Always have a spare CAPIOX FX05 oxygenator for replacement if needed.

- Remove water lines from the old oxygenator module and connect these water lines to the spare oxygenator module. Start water supply to check for leaks.
- Set up and prime spare oxygenator module as described in SET-UP section (page 8).
- When necessary, maintain a decreased patient temperature as prescribed by physician.
- Double-clamp blood inlet and outlet tubing connected to new oxygenator, and cut the tubing between the two clamps. (See Fig.7)
- Stop circulation, double-clamp venous and arterial line connected to old oxygenator, and cut the lines between the two clamps. (See Fig.7)
- Replace old oxygenator with new one by connecting venous and arterial line to the connectors at blood inlet and outlet tubing of new oxygenator. (See Fig.8)

**CAUTIONS**

- After replacement, open recirculation line for debubbling.
- Band and secure all connections in the circuit.

- Start circulation at a low flow rate.
- Connect gas lines removed from old oxygenator module to the new module and start gas supply.
- Connect the purge line of the new oxygenator to the luer port leading to the cardiotomy filter.
- Connect sampling line of new oxygenator to venous side connecting site of old reservoir. (See Fig.9)
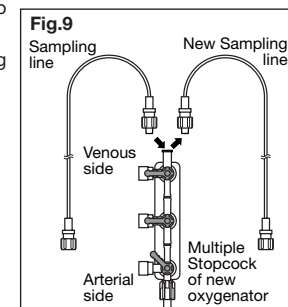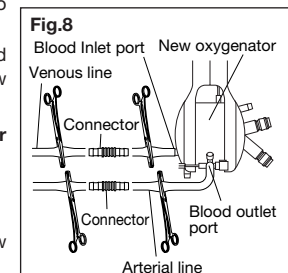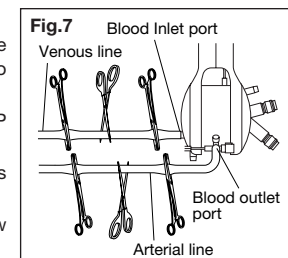
## Vacuum Assisted Venous Drainage (VAVD)

### INTENDED USE

The FX05 Hardshell Reservoir is also intended for use in vacuum assisted venous drainage procedures. The use of a controlled vacuum on the hardshell reservoir improves venous drainage during minimally invasive surgical techniques or regular bypass surgery.

### WARNINGS

- For drug administration in the hardshell reservoir, be aware that negative pressure in this line could draw more drug than intended. To prevent non intended drug amounts to be drawn, it is recommended to include only the intended amounts in the syringe.
- Always allow for a vented reservoir when vacuum is not being applied. Vent by declamping the vacuum release line.
- Avoid obstruction or occlusion of vacuum line to prevent retrograde air to the patient.
- Do not exceed -20 kPa (-150 mmHg) when using vacuum source for assisted venous drainage technique to prevent hemolysis.
- When stopping the pump or when the blood flow rate is low during VAVD, all A-V shunt lines (ex. sampling line, purge line, etc) must be closed to prevent drawing air into the blood side of the oxygenator from the fiber and reversing the blood flow to the reservoir from the arterial side of the patient.
- The roller pump must be occluded properly because it is more likely to draw air into the blood side of the oxygenator from the fiber when VAVD is performed.
- Proceed slowly when returning to atmospheric pressure (elimination of vacuum). A sudden change in pressure may create turbulence of the blood inside the reservoir.
- Do not open the auxiliary port on the venous reservoir during VAVD to prevent bubble formation inside the venous filter.

#### When VAVD is performed in conjunction with a centrifugal pump take caution as below.

- The line between the oxygenator and the centrifugal pump must be clamped before stopping the pump. Not clamping off the arterial line may draw air into the blood side of the oxygenator from the fiber. An arterial line one-way valve is recommended between the oxygenator and the centrifugal pump.
- Caution must be taken because the relationship between blood flow rate and pump speed changes depending on the degree of negative pressure applied.

### PRECAUTIONS


- A moisture trap is required due to the volume of condensation generated.
- A controlled vacuum regulator is required.
- A positive pressure relief valve is required.
- The use of a negative pressure manometer gauge on the hardshell reservoir and a negative pressure relief valve [opening range at -20 kPa (-150 mmHg)] are recommended.
- Please be aware that use of vacuum may cause liquid level to appear higher than actual due to reservoir compliance.

### INSTRUCTIONS FOR USE


#### REQUIRED EQUIPMENT

- |  |   |
|--|---|
| <b>A</b> FX05R venous reservoir  | <b>E</b> Sterile manometer gauge line   |
| <b>B</b> Sterile moisture trap   | <b>F</b> Negative manometer gauge   |
| <b>C</b> Clamp   | <b>G</b> Sterile vacuum release line  |
| <b>D</b> Controlled vacuum regulator [adjustable between 0 to -20 kPa (-150 mmHg)] | <b>H</b> Sterile positive pressure relief valve [opening range from 0 to 1.3kPa (10mmHg)] |
|  | <b>I</b> Sterile gas filter   |

#### SET-UP CONFIGURATION
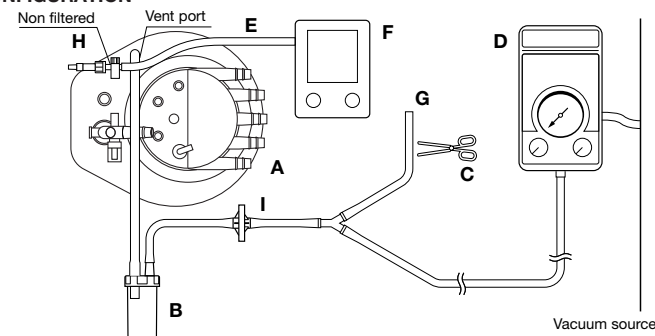
#### STRUCTURE

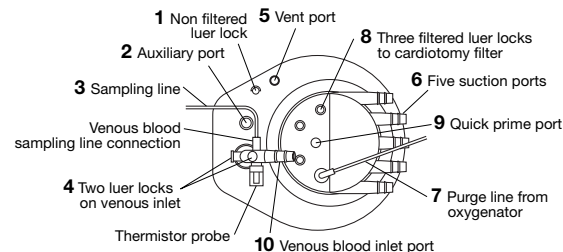Note : 1~10 refers to the Structure diagram with numbered components and A~I refers to the Required Equipment for use with vacuum assisted venous drainage with letter indications.

- Prepare CPB circuit set-up as per standard procedure.
- Verify that blue caps from suction inlet ports **6** and auxiliary port **2** of reservoir are pushed all the way in to insure adequate seal.
- Secure in place all yellow luer caps; all caps are non-vented **4, 8**.
- Install controlled vacuum regulator **D** on wall suction source.
- Attach a 3-way stopcock to the non-filtered luer lock port **1**, and connect the sterile manometer gauge line **E**.
- Connect the negative manometer gauge **F** to the sterile manometer line **E**.
- Also, to the above mentioned (in procedure 5) 3-way stopcock, connect the positive pressure relief valve **H** to the other end of the 3-way stopcock.
- Attach the gas filter **I** to the tubing connected to the moisture trap **B** and the "Y" connector.
- Attach sterile moisture trap **B** tube to the vent port **5** on the hardshell reservoir.

#### <Initiating Bypass>

- Begin with regular venous gravity drainage; at this time, the vacuum release line **G** is unclamped.
- To initiate vacuum assist venous drainage: set the vacuum regulator **D** to -5.3kPa (-40mmHg), then clamp the vacuum release line **G**.
- Monitor the negative pressure inside the hardshell reservoir with the negative manometer gauge **F**.
- Adjust the negative pressure to optimize venous return. Set vacuum regulator **D** between -5.3kPa and -8.0kPa (-40 and -60mmHg).


Note: • The "DO NOT OBSTRUCT" caution label near the vent port is not applicable for vacuum assisted venous drainage.  
• Use an occlusive roller head for the suction and LV vent lines.

#### <Coming off bypass>

Unclamp vacuum release line **G**, venous return will quickly fall. Come off bypass as per standard procedure.

## SAFETY MARGIN

Fig.10 Safety Margin

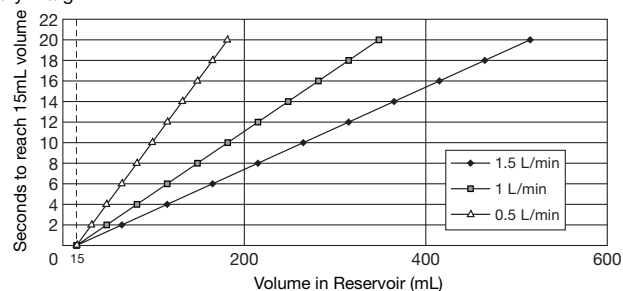\* Minimum operating blood level in CAPiox FX05 Hardshell Reservoir.

The above graph indicates the time in which the blood volume in the reservoir decreases to minimum operating blood level of 15 mL at each flow rate. If the blood volume in reservoir is below 15 mL, bubbles may pass to the oxygenator module.

## PERFORMANCE DATA

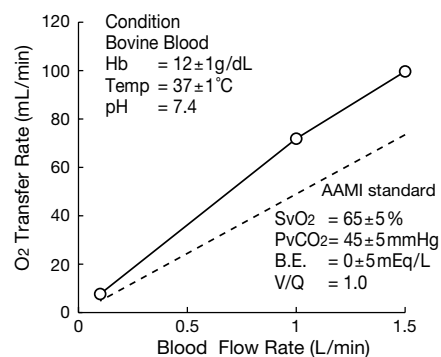
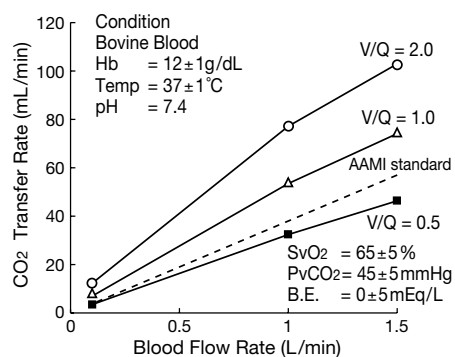
Fig.11 O<sub>2</sub> Transfer Rate (in vitro)Fig.12 CO<sub>2</sub> Transfer Rate (in vitro)

Fig.13 Heat Exchanger Performance Factor (in vitro)


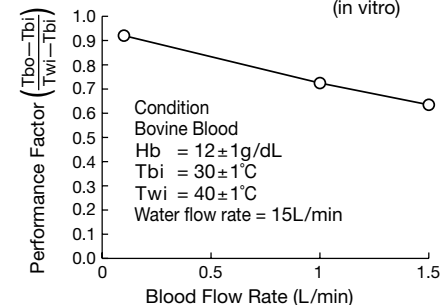
Fig.14 Blood Side Pressure Drop (in vitro)
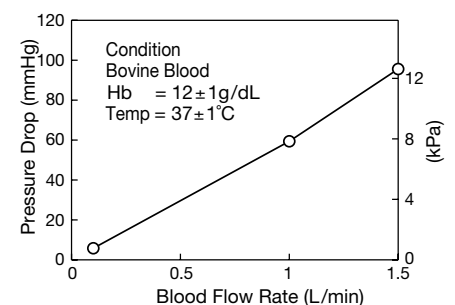
## SPECIFICATIONS

## Product Description

Reorder Number	Description
CX*FX05RW	Oxygenator with integrated arterial filter with Hardshell Reservoir
CX*FX05RE	Oxygenator with integrated arterial filter with Hardshell Reservoir
CX*FX05W	Oxygenator with integrated arterial filter
CX*FX05E	Oxygenator with integrated arterial filter

(Refer to "Direction of blood outlet port", page 3)

## CAPiox FX05 Oxygenator Module

COMPONENT	SPECIFICATIONS	
Housing	Material	Polycarbonate
Fibers	Material Surface area	Microporous Polypropylene Approx. 0.5 m <sup>2</sup>
Arterial filter	Material	Polyester screen type Pore size 32 μm
Heat Exchanger	Material Surface area	Stainless steel Approx. 0.035 m <sup>2</sup>
Blood Flow Range		Min. 0.1 L/min Max. 1.5 L/min
Reference Blood Flow (AAMI standard)		2.5 L/min
Priming Volume ( Static )		43 mL
Blood inlet port ( from pump )		1/4" (6.4mm)
Blood outlet port		1/4" (6.4mm)
Luer port (for recirculation or blood cardioplegia)		one luer lock on blood outlet port
Gas inlet port		1/4" (6.4mm)
Gas outlet port		5/16" (7.9mm)
Water ports		1/2" (12.7mm) Hansen quick connect fittings
Maximum Pressure		Blood inlet 133kPa (1,000 mmHg) Water inlet 196kPa (2 kgf/cm <sup>2</sup> ) (28.5 psi)

## Hardshell Reservoir Section

COMPONENT	SPECIFICATIONS	
Housing	Material	Polycarbonate
Blood Flow Range		Min. 0.1 L/min Max. 1.5 L/min
Blood Storage Capacity		1,000 mL
Minimum Operating Volume		15 mL
Venous Filter	Material	Polyester screen type Pore size 47 μm
Cardiotomy Filter	Material	Polyester depth type
Defoaming Part	Material	Polyurethane foam
Venous blood inlet port		1/4" (6.4mm) rotatable
Blood outlet port ( to pump )		1/4" (6.4mm)
Suction ports		five 3/16"—1/4" (4.8mm—6.4mm) rotatable
Quick prime port		1/4" (6.4mm)
Vent port		1/4" (6.4mm)
Auxiliary port		1/4"—3/8" (6.4mm—9.5mm)
Luer ports		• three filtered luer locks to cardiotomy filter • a non-filtered luer lock • two luer locks on venous inlet
Maximum Sustainable Negative Pressure in Reservoir		-20 kPa (-150 mmHg)

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## LIMITED WARRANTY

TERUMO Corporation warrants that the product has been manufactured and packaged with reasonable care, and will be free from defects under normal and proper use. TERUMO Corporation will not be liable for any direct, indirect, incidental, special, general or consequential loss, damage, or expense from the use of this product. TERUMO Corporation's sole obligation under this warranty shall be to repair or replace, at its option, any products that it believes were defective at the time of shipment if notice thereof is received before the expiration date indicated on the package. THIS WARRANTY IS THE SOLE AND ONLY WARRANTY OF TERUMO CORPORATION AND IS EXPRESSLY IN LIEU OF AND EXCLUDES ANY OTHER EXPRESSED OR IMPLIED WARRANTY OF ANY KIND, INCLUDING THE WARRANTY OF MERCHANTABILITY OR FITNESS FOR ANY PARTICULAR PURPOSE. IN ADDITION, THE REMEDIES SET FORTH HEREIN ARE THE SOLE REMEDIES AVAILABLE TO ANY PERSON OR ENTITY. No agent, employee or representative of TERUMO Corporation has any authority to change, amend or alter the foregoing or assume or bind TERUMO Corporation to any additional liability or responsibility in connection with this product.

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**TERUMO EUROPE N.V.**  
INTERLEUVENLAAN 40, 3001 LEUVEN, BELGIUM

**TERUMO CORPORATION**  
44-1, 2-CHOME, HATAGAYA, SHIBUYA-KU, TOKYO 151-0072, JAPAN  
MADE IN JAPAN

® : Registered Trademark	BABY-FX and Xcoating are trademarks of TERUMO CORPORATION.	©TERUMO CORPORATION
® : Marque Enregistrée	BABY-FX et Xcoating sont des marques de TERUMO CORPORATION.	Jan. 2009 09A15
® : Eingetragenes Warenzeichen	BABY-FX und Xcoating sind Marken der TERUMO CORPORATION.	
® : Marca Registrada	BABY-FX y Xcoating son marcas de TERUMO CORPORATION.	
® : Marchio Registrato	BABY-FX e Xcoating sono marchi della TERUMO CORPORATION.	
® : Geregistreerd Handelsmerk	BABY-FX en Xcoating zijn handelsmerken van TERUMO CORPORATION.	

**CAPIOX FX05**  
**CX\_FX05\_M6\_50\_003**

CF:10036 (pink)

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ANN ARBOR, MICHIGAN 48103-9300, U.S.A.