





Instruction Manual

CO₂-1 SMA

Specification:

Fiber optic CO₂ meter for use with CO₂ sensor spots, flow-through cells & dipping probes

Software:

PreSens Measurement Studio (PMS2)

Document filename: IM_CO2-1 SMA_dv1

All rights reserved. No parts of this work may be reproduced in any form or by any means - graphic, electronic, or mechanical, including photocopying, recording, taping, or information storage and retrieval systems - without the written permission of the publisher.

Products that are referred to in this document may be either trademarks and/or registered trademarks of the respective owners. The publisher and the author make no claim to these trademarks.

While every precaution has been taken in the preparation of this document, the publisher and the author assume no responsibility for errors or omissions, or for damages resulting from the use of information contained in this document or from the use of programs and source code that may accompany it. In no event shall the publisher and the author be liable for any loss of profit or any other commercial damage caused or alleged to have been caused directly or indirectly by this document.

Specifications may change without prior notice.

Manufacturer

PreSensAm BioPark 11Phone +49 941 94272100info@PreSens.dePrecision Sensing GmbH93053 Regensburg, GermanyFax +49 941 94272111www.PreSens.de

Table of Contents

1	Pretace	
2	Safety Notes	8
3	Description of the CO2-1SMA	9
4	Set-Up	10
4.1	Connecting the Devices	10
4.2	Important Considerations for USB Handling	12
5	Technical Data	13
6	Operational Notes	14
6.1	Optical Output	
6.2	Maintenance	14
6.3	Service	14
7	Concluding Remarks	15

1 Preface

You have chosen a new, innovative technology for measuring carbon dioxide.

The CO_2 -1 SMA is a compact fiber optic CO_2 meter. It is based on a novel technology, which creates very stable, internally referenced measured values. This allows a more flexible use of CO_2 sensors in various fields of interest.

Optical CO₂ sensors (also called optrodes) have several important features:

- They are small.
- Their signal does not depend on the flow rate of the sample.
- They can be physically divided from the measuring system which allows a noninvasive measurement.
- They can be used in disposables.

Therefore, they are ideally suited for the examination of small sample volumes, for highly parallelized measurements in disposables, and for biotechnological applications. A set of different CO₂ minisensors, flow-through cells and non-invasive sensors is available to make sure you have the sensor which matches your application.

Please feel free to contact our service team to find the best solution for your application.

Your PreSens Team

PLEASE READ THE FOLLOWING INSTRUCTIONS CAREFULLY BEFORE WORKING WITH THIS DEVICE. WHEN DISREGARDING THESE INSTRUCTIONS THE SAFETY OF THE DEVICE CAN BE IMPAIRED.

2 Safety Notes

It is the customer's responsibility to validate the sensor and transmitter under end-user conditions according to safety precautions of the application to ensure that the use of the sensor is safe and suitable for the intended purpose.

PreSens is explicitly not liable for direct or indirect losses caused by the application of these measurement systems. In particular it has to be considered that malfunctions can occur due to the naturally limited lifetime of the sensor depending on the respective application. The set-up of backup measurement stations is recommended when using the sensors in critical applications to avoid consequential losses. It is the customer's responsibility to install a suitable safety system in the event of sensor failure.

3 Description of the CO₂-1SMA



Fig. 1 CO₂-1 SMA optical carbon dioxide meter

The single channel CO₂-1 SMA device is compatible with non-invasive sensors, dipping probes and flow-through cells of type CD1 (1 - 25 % CO₂). The device is USB powered and does not need an extra power adapter. CO₂-1 SMA has temperature compensation, so even in environments with changing temperatures precise CO₂ measurements can be performed. The CO₂ meter is operated with the PreSens Measurement Studio 2 software, which enables simultaneous control of several devices. With numerous features and additional pressure compensation, the software makes the CO₂-1 SMA suitable for numerous applications.

Several single- and multichannel devices can be connected to the PreSens Measurement Studio 2 software, so you can measure with as many sensors simultaneously as required.

4 Set-Up

4.1 Connecting the Devices

1. Remove the protective cap from the male plug on the polymer optical fiber (or the probe) and insert it in the SMA connector of the CO₂-1 SMA. The safety nut must be screwed on.



Fig. 1 Connecting the SMA plug to the CO₂-1 SMA

For further information on sensor handling and calibration please refer to the respective sensor instruction manual.

Connect the USB cable to the connector on the CO₂-1 SMA and to a USB port of your PC / notebook.



Fig. 2 Connecting the USB cable

3. If required, you can now connect a Pt100 temperature sensor to the device.



Fig. 3 Connecting the temperature sensor to CO₂-1 SMA

There is a red mark on the temperature sensor connector of the CO₂-1 SMA as well as on the temperature sensor plug. Match those two marks before inserting the temperature sensor plug into the connector on the CO₂-1 SMA; else, the plug might be damaged.

4. Now start the PreSens Measurement Studio (PMS2) software on your PC / notebook. (The CO₂-1 SMA always has to be connected via USB before the software is started.) After successful initialization, the software main screen is displayed and the CO₂-1 SMA will show in the **Device** section.

For further details on how to use the PreSens Measurement Studio please refer to the software instruction manual.

4.2 Important Considerations for USB Handling

- PreSens recommends the use of a dedicated USB 2.0 PCI Card to connect and handle USB PreSens Devices with a Desktop PC.
- In order to enhance the system stability avoid the use of USB hubs and connect PreSens Devices directly to your PC / notebook USB Ports.
- If possible, disconnect all other USB devices that are not in use, as they may reduce or disturb the USB resources of your PC / notebook.
- Docking stations may also reduce or disturb the USB resources of your PC / notebook and therefore affect the correct function of the software.
- It is also recommended to disable the Power Saving Settings of your USB Root controller.

5 Technical Data

OPTICAL SENSOR	
CO ₂ sensor	CD1
Optical connector	SMA
Channels	1; up to 10 CO ₂ -1 SMA can be controlled simultaneously via PMS2 software
LED peak wavelength	470 nm

TEMPERATURE SENSOR		
Potentiometric temperature sensor	Range	0 – 50 °C
(Pt100)	Resolution	± 0.1 °C

POWER SUPPLY	
Supply voltage	5 VDC (USB-2.0-Mini-B)
Current / Power	400 mA

DIGITAL INTERFACE	
USB interface cable to PC	Cable included

ENVIRONMENTAL CONDITIONS	
Operating temperature	0 °C to 50 °C
Storage temperature	- 20 °C to 70 °C
Relative humidity	0 % to 80 % (non-condensing)

DIMENSIONS / WEIGHT	
HxLxW	ca. 30 mm x 95 mm (with connectors) x 34 mm
Weight	128 g

6 Operational Notes

6.1 Optical Output

The SMA connector is a high precision optical component. Please keep it clean and dry. Always use the rubber cap to close the output when not in use.

6.2 Maintenance

The transmitter is maintenance-free.

The housing should be cleaned with a cloth only. Avoid any moisture entering the housing. Never use benzine, acetone, alcohol or any other organic solvents.

The SMA fiber connector of the sensor can be cleaned with lint-free cloth or a cleaning implement for SMA connectors only.

6.3 Service

Alignment, rework or repair work may only be carried out by the manufacturer:

PreSens Precision Sensing GmbH

Am BioPark 11 93053 Regensburg Germany

Phone +49 941 94272100 Fax +49 941 94272111

info@PreSens.de www.PreSens.de

Please contact our service team in case of any question. We look forward to helping you and are open for any proposition or criticism.

7 Concluding Remarks

Dear Customer,

With this manual, we hope to provide you with an introduction to work with the CO_2 -1 SMA fiber optic CO_2 meter.


This manual does not claim to be complete. We are endeavored to improve and supplement this version.

We are looking forward to your critical review and to any suggestions you may have.

You can find the latest version at www.PreSens.de.

With best regards,

Your PreSens Team

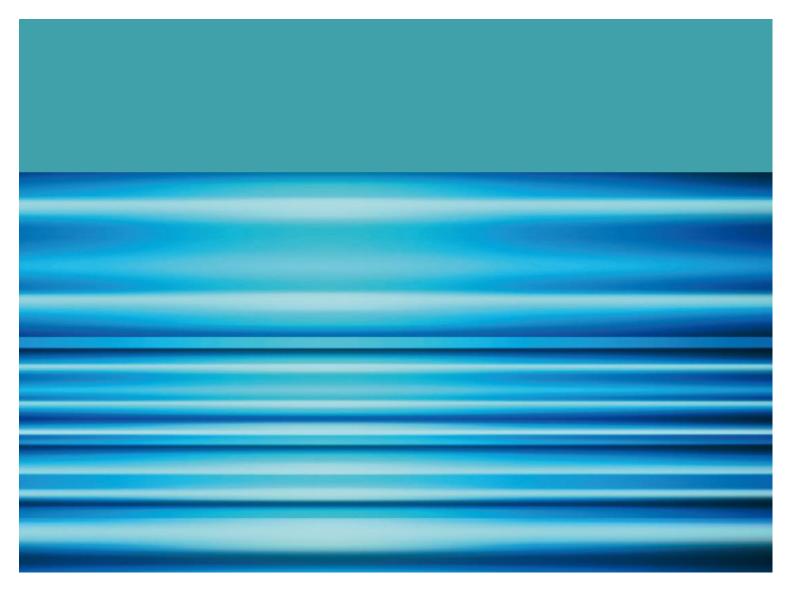

Manufacturer

PreSens Precision Sensing GmbH

Am BioPark 11 93053 Regensburg Germany

Phone +49 941 94272100 Fax +49 941 94272111

info@PreSens.de www.PreSens.de