## pH-1 SMA HP5 / LG1









## pH-1 SMA HP5 / LG1

Specification:

Fiber optic pH meters for use with pH sensor spots, flow-through cells & dipping probes

Software:

PreSens Measurement Studio (PMS2)

Document filename: IM\_pH-1 SMA\_dv1

All rights reserved. No parts of this work may be reproduced in any form or by any means - graphic, electronic, or mechanical, including photocopying, recording, taping, or information storage and retrieval systems - without the written permission of the publisher.

Products that are referred to in this document may be either trademarks and/or registered trademarks of the respective owners. The publisher and the author make no claim to these trademarks.

While every precaution has been taken in the preparation of this document, the publisher and the author assume no responsibility for errors or omissions, or for damages resulting from the use of information contained in this document or from the use of programs and source code that may accompany it. In no event shall the publisher and the author be liable for any loss of profit or any other commercial damage caused or alleged to have been caused directly or indirectly by this document.

Specifications may change without prior notice.

Manufacturer

PreSensAm BioPark 11Phone +49 941 94272100info@PreSens.dePrecision Sensing GmbH93053 Regensburg, GermanyFax +49 941 94272111www.PreSens.de

## **Table of Contents**

1	Preface	7
2	Safety Notes	8
3	Description of the pH-1 SMA HP5 and pH-1 SMA LG1	9
4	Set-Up	10
5	Technical Data	12
6 6.2 6.3	Operational Notes Maintenance Service	13
7	Concluding Remarks	14

#### 1 Preface

You have chosen a new, innovative technology for measuring pH.

The pH-1 SMA is a compact fiber optic pH meter. It is based on a novel technology, which creates very stable, internally referenced measured values. This allows a more flexible use of pH sensors in various fields of interest.

Optical pH sensors (also called optrodes) have several important features:

- They are small.
- Their signal does not depend on the flow rate of the sample.
- They can be physically divided from the measuring system which allows a noninvasive measurement.
- They can be used in disposables.

Therefore, they are ideally suited for the examination of small sample volumes, for highly parallelized measurements in disposables, and for biotechnological applications. A set of different pH minisensors, flow-through cells and non-invasive sensors is available to make sure you have the sensor which matches your application.

Please feel free to contact our service team to find the best solution for your application.

Your PreSens Team

PLEASE READ THE FOLLOWING INSTRUCTIONS CAREFULLY BEFORE WORKING WITH THIS DEVICE. WHEN DISREGARDING THESE INSTRUCTIONS THE SAFETY OF THE DEVICE CAN BE IMPAIRED.

#### 2 Safety Notes

It is the customer's responsibility to validate the sensor and transmitter under end-user conditions according to safety precautions of the application to ensure that the use of the sensor is safe and suitable for the intended purpose.

PreSens is explicitly not liable for direct or indirect losses caused by the application of these measurement systems. In particular it has to be considered that malfunctions can occur due to the naturally limited lifetime of the sensor depending on the respective application. The set-up of backup measurement stations is recommended when using the sensors in critical applications to avoid consequential losses. It is the customer's responsibility to install a suitable safety system in the event of sensor failure.

# 3 Description of the pH-1 SMA HP5 and pH-1 SMA LG1



Fig. 1 pH-1 SMA HP5 and pH-1 SMA LG1

Due to their small size these fiber optic pH meters can be set up almost anywhere. As their names already indicate pH-1 SMA HP5 is compatible with all sensors type HP5 / HP8 (measurement range pH 5.5-8.5) and pH-1 SMA LG1 is compatible with LG1 sensors (measurement range pH 4.0-7.5). Both pH meters have temperature compensation, so even in environments with changing temperature precise pH measurements can be performed. They are USB-powered so there is no need for extra power connection, cables and adapters. The pH-1 SMA devices are controlled with the PreSens Measurement Studio 2 software (free download on <a href="https://www.presens.de/support-services/download-center/software.html">https://www.presens.de/support-services/download-center/software.html</a> ), which enables simultaneous control of several PreSens pH,  $O_2$  and  $CO_2$  devices. With numerous features, the software makes pH-1 SMA HP5 / LG1 the ideal pH measurement tool for almost any application.

#### 4 Set-Up

1. Remove the protective cap from the male plug on the polymer optical fiber (or the probe) and insert it in the SMA connector of the pH-1 SMA. The safety nut has to be screwed on.



Fig. 1 Connecting the SMA plug to the pH-1 SMA

For further information on sensor handling and calibration please refer to the respective sensor instruction manual.

Connect the USB cable to the connector on the pH-1 SMA and to an USB port of your PC / notebook.



Fig. 2 Connecting the USB cable

3. If required, you can now connect a Pt100 temperature sensor to the device (not supplied, must be ordered separately). There is a red mark on the temperature sensor connector of the pH-1 SMA as well as on the temperature sensor plug. Match those two marks before inserting the temperature sensor plug into the connector on the device, else the plug might get damaged.



Fig. 3 Connecting the temperature sensor

4. Now start the PreSens Measurement Studio (PMS2) software on your PC / notebook. (The pH-1 SMA always has to be connected via USB before the software is started.) After successful initialization the software main screen is displayed and the pH-1 SMA will show in the **Device** section.

For further details on how to use the PreSens Measurement Studio please refer to the software instruction manual (<a href="https://www.presens.de/support-services/download-center/manuals-getting-started.html">https://www.presens.de/support-services/download-center/manuals-getting-started.html</a>).

#### Important Considerations for USB Handling

- PreSens recommends the use of a dedicated USB 2.0 PCI Card to connect and handle USB PreSens Devices with a Desktop PC.
- In order to enhance the system stability avoid the use of USB hubs and connect PreSens Devices directly to your PC / notebook USB Ports.
- If possible, disconnect all other USB devices that are not in use, as they may reduce or disturb the USB resources of your PC / notebook.
- Docking stations may also reduce or disturb the USB resources of your PC / notebook and therefore affect the correct function of the software.
- It is also recommended to disable the Power Saving Settings of your USB Root controller.

### **5 Technical Data**

OPTICAL SENSOR	
pH sensor	pH-1 SMA HP5: HP5 / HP8
	pH-1 SMA LG1: LG1
Optical connector	SMA compatible, 2 mm PMMA fiber
Channels	1; up to 10 pH-1 SMA can be controlled simultaneously
Chameis	via PMS2 software
LED peak wavelength	pH-1 SMA HP5: 470 nm
LLD peak wavelengin	pH-1 SMA LG1: 505 nm

TEMPERATURE SENSOR		
Potentiometric temperature sensor	Range	0 °C to + 50 °C
(Pt100)	Resolution	± 0.1 °C

POWER SUPPLY	
Supply voltage	5 VDC (USB-2.0-Mini-B)
Current / Power (peak value)	400 mA

DIGITAL INTERFACE	
USB interface cable to PC	Cable included

ENVIRONMENTAL CONDITIONS	
Operating temperature	0 °C to + 50 °C
Storage temperature	- 20 °C to + 70 °C
Relative humidity	0 % to 80 % (non-condensing)

DIMENSIONS / WEIGHT	
HxLxW	ca. 30 x 95 (with connectors) x 34 mm
Weight	128 g

#### **6 Operational Notes**

The SMA connector is a high precision optical component. Please keep it clean and dry. Always use the rubber cap to close the output when not in use.

#### 6.2 Maintenance

The transmitter is maintenance-free.

The housing should be cleaned with a cloth only. Avoid any moisture entering the housing. Never use benzine, acetone, alcohol or any other organic solvents.

The SMA fiber connector of the sensor can be cleaned with lint-free cloth or a cleaning implement for SMA connectors only.

#### 6.3 Service

Alignment, rework or repair work may only be carried out by the manufacturer:

## PreSens Precision Sensing GmbH

Am BioPark 11 93053 Regensburg Germany

Phone +49 941 94272100 Fax +49 941 94272111

info@PreSens.de www.PreSens.de

Please contact our service team in case of any question. We look forward to helping you and are open for any proposition or criticism.

### 7 Concluding Remarks

Dear Customer,

With this manual, we hope to provide you with an introduction to work with the pH-1 SMA HP5 / LG1 fiber optic pH meters.


This manual does not claim to be complete. We are endeavored to improve and supplement this version.

We are looking forward to your critical review and to any suggestions you may have.

You can find the latest version at www.PreSens.de.

With best regards,

Your PreSens Team

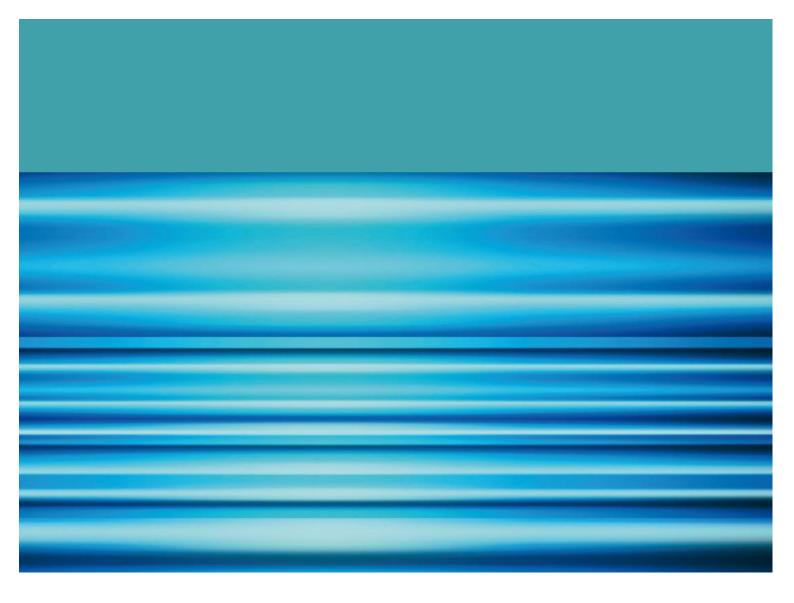

#### Manufacturer

PreSens Precision Sensing GmbH

Am BioPark 11 93053 Regensburg Germany

Phone +49 941 94272100 Fax +49 941 94272111

info@PreSens.de www.PreSens.de